




The Talon® Splint from Space Maintainers

- Will not leach out plasticizers into the oral cavity.
- Based on clear, thermo-elastic material with memory.
- Combines fit and comfort of a soft splint with the wear resistance, accuracy and adjustability of a hard acrylic occlusal surface.
- Base becomes flexible when placed under hot tap water for easy insertion.
- As it cools, splint returns to original shape, conforming to shape of patient's teeth (allows for use of natural buccal undercuts and lingual interproximal spaces for reliable retention).
- Clasp for retention almost eliminated.
- Can be utilized for all splint designs (e.g., Gelb, MORA, May, Sears, Jankelson, Tanner, etc.)

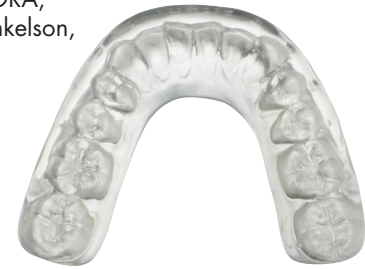

Visit us at www.aurumgroup.com

Local Pick-up and Delivery Service Available In:

SPACE MAINTAINERS LABORATORIES CANADA LTD.

CALGARY

115 - 17TH AVENUE S.W., CALGARY, AB T2S 0A1
(403) 228-5120 • TOLL FREE 1-800-661-1169
spacemaint@aurumgroup.com

EDMONTON

11007 - 106TH AVENUE, EDMONTON, AB T5H 4R7
(780) 423-1904 • TOLL FREE 1-800-661-2745
aurumedm@aurumgroup.com

SASKATOON

336 - 6TH AVENUE NORTH, SASKATOON, SK S7K 2S5
(306) 665-8815 • TOLL FREE 1-800-665-8815
aurumsask@aurumgroup.com

VANCOUVER

936 WEST 8TH AVENUE, VANCOUVER, BC V5Z 1E5
(604) 737-2010 • TOLL FREE 1-800-663-1721
aurumvan@aurumgroup.com

VICTORIA

1928 OAK BAY AVENUE, VICTORIA, BC V8R 1C9
(250) 595-2314 • TOLL FREE 1-800-663-6364
aurumvic@aurumgroup.com

OTTAWA

1175 CECIL AVENUE, OTTAWA, ON K1H 7Z6
(613) 736-1946 • TOLL FREE 1-800-267-7040
classcottawa@aurumgroup.com

TORONTO

40 PIPPIN ROAD, UNIT 11 & 12, CONCORD, ON L4K 4M6
(416) 410-1330 • TOLL FREE 1-800-268-4294
classictoronto@aurumgroup.com



A Canadian Group of Companies

© Space Maintainers Laboratories Canada (2010).
Reproduction of this work in whole or in part, by any means whatsoever,
is strictly prohibited without the express written consent of
Space Maintainers Laboratories Canada. All rights reserved.

The Talon® Splint

Safest Splint in Laboratory Testing!



The Talon[®] Splint

*Proven Performance...
No Discernable Leaching
of Free Monomer!*



Outstanding Performance Proven in Laboratory Testing

- Will not leach out plasticizers into the oral cavity.
- Gas chromatographic studies verify residual monomer level of only 0.14%.
- Will not harden over time, has "memory" and will not warp under normal wear and handling of the appliance.

Comfortable Retentive Base for Trouble-free, Instant Fit

- Created from soft, thermoplastic, resilient Talon[®] polymer.
- Becomes flexible when heated under warm tap water.

- Returns to original shape as it cools, utilizing natural buccal undercuts and lingual interproximal spaces for retention.
- Eliminates time spent trying to find "tight spots" on appliance due to pressure points.
- Maintains its original flexibility for years (unlike other thermoplastics currently available).

Accuracy and Adjustability of a Hard Occlusal Surface

- Hard acrylic chemically bonded to processed Talon[®] material to form splint's occlusal surface. Allows easy, precise occlusal adjustments.
- Proper proprioceptive response with the opposing dentition to be maintained.

Wide Range of Applications

- Great treatment option for bruxism, periodontal splinting, orthodontic stabilization, or TMJ dysfunction.
- Can be utilized for all splint designs (i.e., Gelb, MORA, May, Sears, Jankelson, Tanner, etc.)

For splints that fit and function exactly as you prescribed, without leaching monomer, talk to the professionals at your closest Space Maintainers laboratory (please see reverse for locations and TOLL FREE contact information).

SPACE MAINTAINERS[®]
LABORATORIES CANADA

Visit us at
www.aurumgroup.com